


UroCure


ArcTV™

TRANSVAGINAL SLING SYSTEM

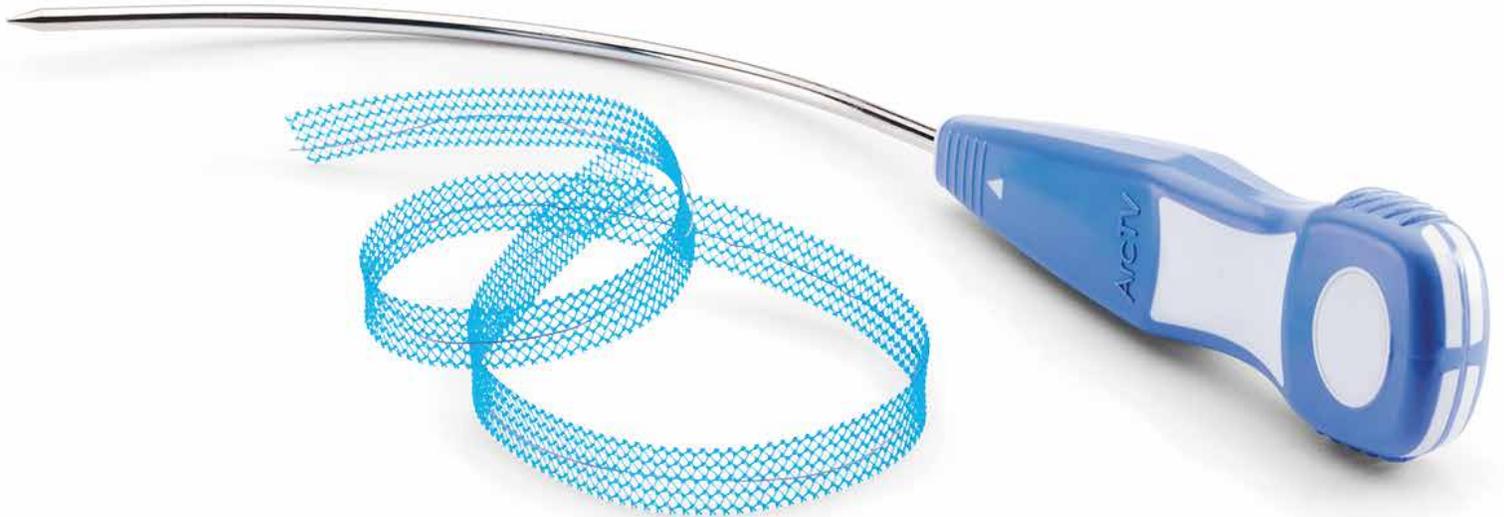

INTENTIONALLY DESIGNED. THOUGHTFULLY REFINED.

Developed in collaboration with leading surgeons, ArcTV is designed and engineered for the safe and accurate placement of a mid-urethral sling by the transvaginal route.

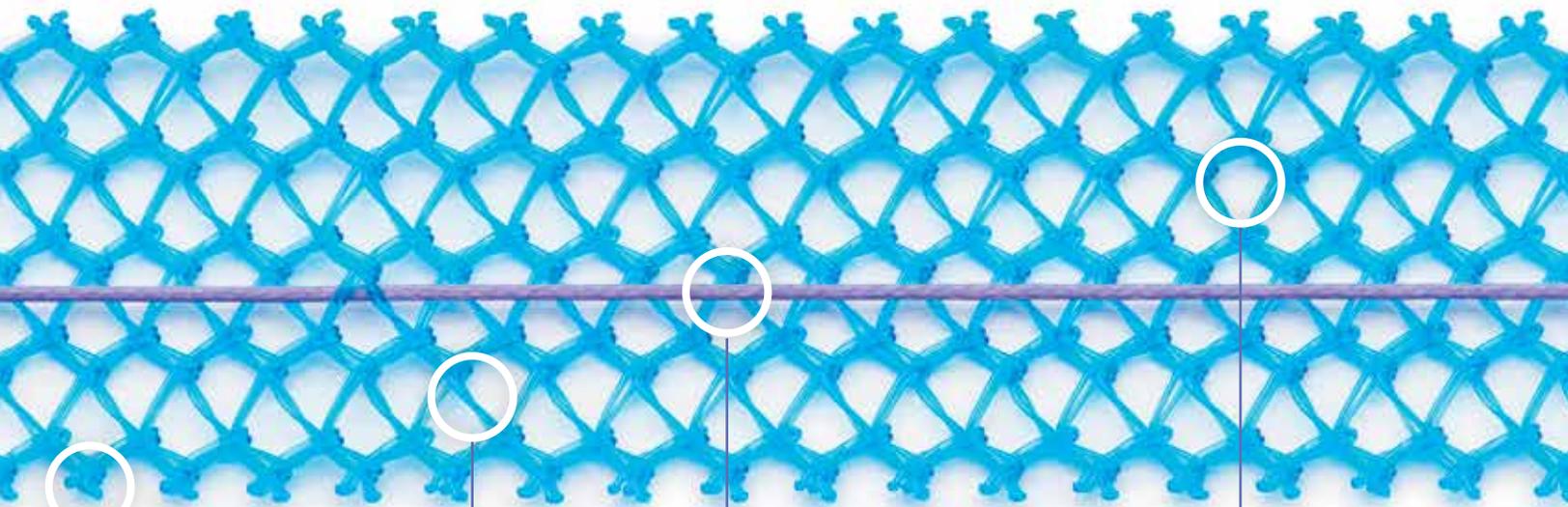

**EQUI-CUT™ SLING**

has a uniform and smooth edge that is designed to reduce tissue trauma during sling integration.

**BLUE SLING**

is designed to improve visibility for cystoscopy and any future intervention.


**PATENTED ABSORBABLE STABILIZING SUTURE**

uniquely maintains the sling's architecture at the midurethra and is designed to resist deformation during sling placement and sheath removal.


**KNITTED SLING PATTERN**

is designed with large pore sizes to promote tissue in-growth.<sup>1</sup>

# PRODUCT OVERVIEW

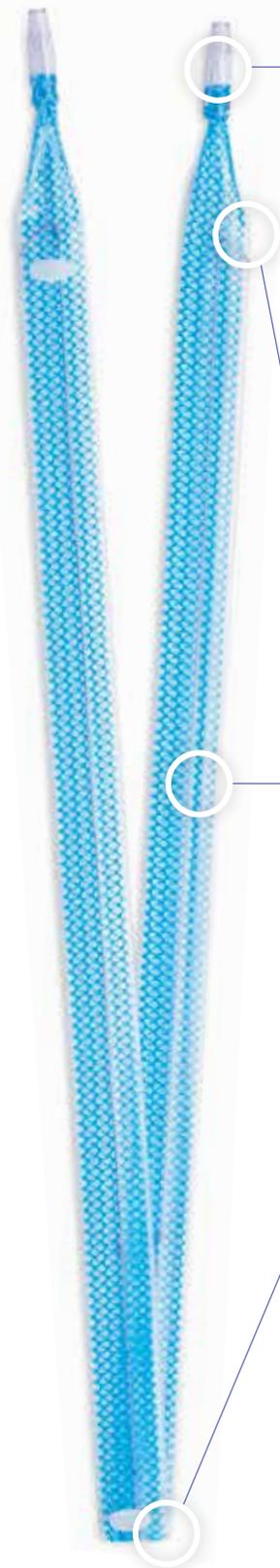


○ **"CLEAR-CLICK" RAPID CONNECTOR** simplifies easy and secure sling connection to the needles.

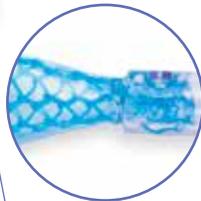

○ **ATRAUMATIC TAPERED SHEATHS** are designed to minimize tissue trauma and decrease the force required during sheath-sling passage.

○ **EASILY-REMOVABLE SHEATHS** designed to maintain sling integrity during removal.

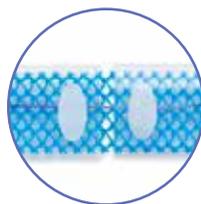

○ **MIDLINE DOTS AND NON-OVERLAPPING SHEATHS** indicate the center of the sling for accurate placement and help with recognizing sling twisting.



○ **QUICK-RELEASE BUTTON** disengages the handle to allow reapplication for the second needle passage.



○ **ERGONOMICALLY-DESIGNED HANDLE & NEEDLE** are specifically balanced to optimize needle control throughout the entire route of passage.

○ **CAREFULLY-DESIGNED 4mm-DIAMETER NEEDLES** offer appropriate rigidity with minimal needle deflection. The 110° arc is designed to follow the trajectory along the back of the pubic bone.

○ **CONICAL NEEDLE TIP** is designed to provide smooth needle passage to minimize tissue trauma.


○ **ALIGNMENT MARKINGS** provide a visual reference during needle passage.

# ORDERING INFORMATION:

ORDERING PART NUMBER:	PRODUCT:	QUANTITY:
A-TV-1001	UROCURE ArcTV TRANSVAGINAL SLING SYSTEM	INCLUDES ONE (1) ArcTV SLING SYSTEM

Note: Product shipped in quantity of two to support two sling procedures.

If you would like more information, please call Customer Service at (612) 466-0117 or email us at [customerservice@urocure.com](mailto:customerservice@urocure.com).

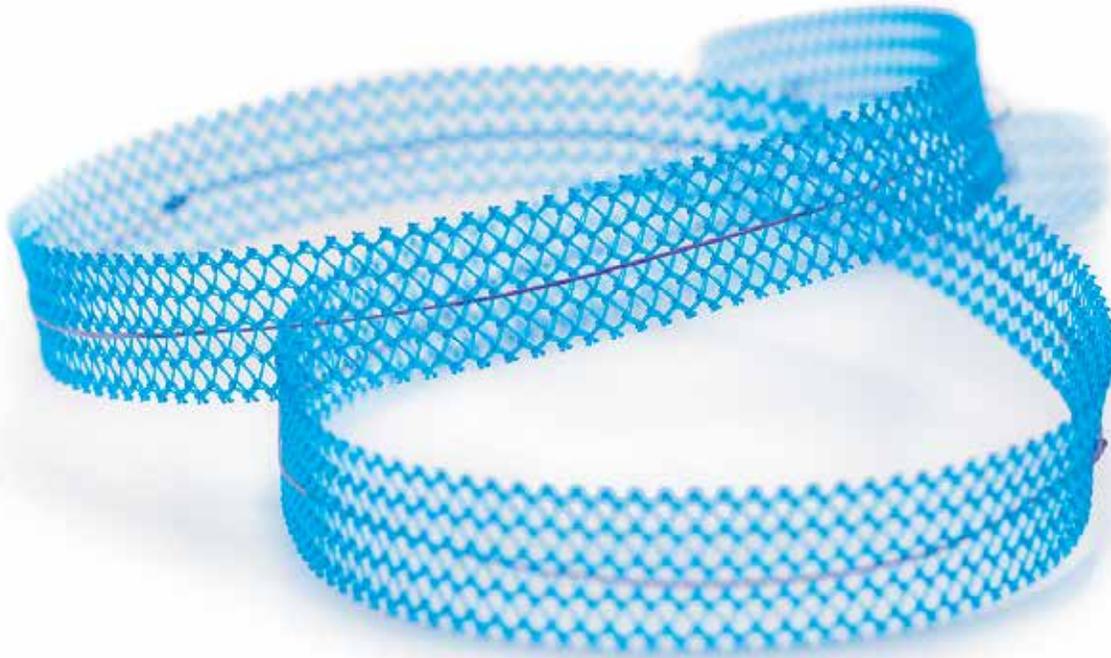

**CAUTION:** U.S. federal law restricts this device to sale by or on the order of a physician.

**CAUTION:** This device is intended for use only by physicians with adequate training and experience in the surgical treatment of stress urinary incontinence. Indications, contraindications, warnings and precautions, adverse events and instructions for use can be found in the ArcTV Instructions for Use.



# A NEW WAY FORWARD FOR WOMEN'S HEALTH

UroCure is a new company formed to return to market a proven medical device and procedure for treating stress urinary incontinence in women. UroCure's ArcTV Sling incorporates popular features of the AMS vaginal retropubic sling—particularly the handle and needle assembly—with enhancements such as a blue sling for improved visibility and softened laser-cut edges. Every aspect of the ArcTV has been carefully considered for how it serves both doctor and patient.

Led by John Nealon, President and CEO, and Dr. David Staskin, Chief Medical Officer, UroCure's leadership team has over 100 years of collective experience in the field of women's health. Dr. Staskin's patented designs for midurethral slings helped position AMS as a global leader in urogynecological products, starting in the early 2000s. Dr. Staskin and Mr. Nealon—together with other medical and industry experts—collaborated to bring an enhanced version of Dr. Staskin's patented design to the marketplace to help women needing surgical intervention for stress urinary incontinence.

As UroCure launches the ArcTV, our company is taking a three-pronged approach to promote transparency and high quality outcomes for your patients. UroCure is working with the American Urogynecological Society (AUGS) and their AQUIRE registry to follow patients post-surgery. In addition to the registry, we have impaneled a Quality and Safety Oversight Committee (QSOC) of leading academic physicians that will monitor those outcomes and provide guidance to UroCure's activities moving forward. And finally, we have created a product training protocol to be sure physicians have the knowledge necessary to successfully implant the ArcTV sling.

UroCure sees itself in partnership with physicians and their patients. We have intentionally incorporated feedback in every step of our process—from the design of the ArcTV, its introduction to hospitals and physicians, and to the surgical procedure. This clear and transparent development process and market approach is groundbreaking and demonstrates a new way forward to maximize patient outcomes.

# UroCure

ADVANCING WOMEN'S HEALTH

© 2019 UroCure L.L.C. All rights reserved. UroCure, ArcTV, Equi-Cut and its other trademarks, logos and slogans are the intellectual property of UroCure, L.L.C. and may not be used without permission. | LBL-1001-001 Rev A | April 2019

701 N 3rd St, Suite 110 | Minneapolis, MN 55401  
(612) 466-0117 | [customerservice@urocure.com](mailto:customerservice@urocure.com) | [www.urocure.com](http://www.urocure.com)